


**THE MEDSTAR™ 300  
MONITORING SYSTEM**



**Simple monitoring system  
for blood pressure, weight, glucose,  
and a fourth parameter or device**

**CYBERNET MEDICAL®**

A Revolution In Home Monitoring For Chronic Disease Patients

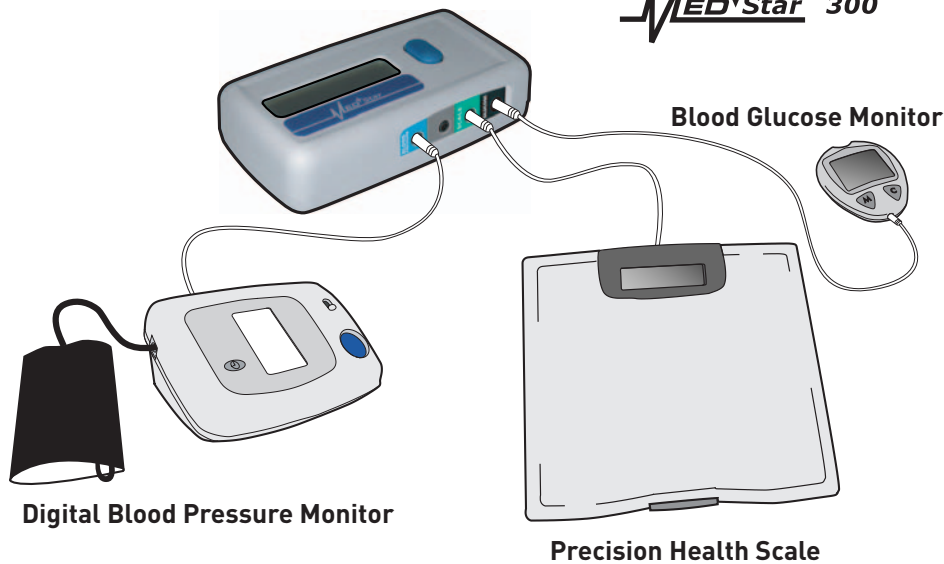

## MEDSTAR SPECIFICATIONS:

Size: 4.94" x 2.75" x 1.28"

Weight: 5.7 Ounces with Batteries

Battery: 2 AA Alkalines


Battery Life: 6 Months

## LCD DISPLAY:

Text: 2 Lines x 16 Characters Alphanumeric

## PHONE COMMUNICATIONS TO MEDSTAR SERVER

## OPTIONAL WEB-BASED DATA MANAGEMENT SOFTWARE

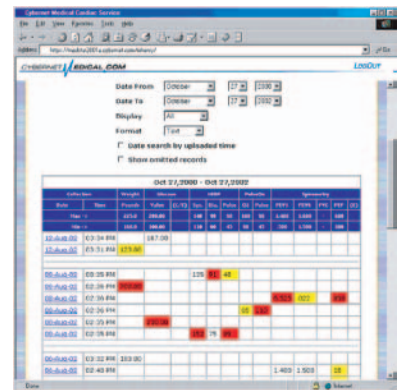

- Simple patient operation
- Electronic serial number for patient identification
- Durable, water resistant, and small in physical size
- Independent battery power
- Stores data from device reading until successfully transferred
- Records time and date of all measurements

The MedStar 300 is currently optimized for use with the following measurement devices. Other devices may be supported. Please see [www.cybernetmedical.com](http://www.cybernetmedical.com) for details.\*

	Manufacturer	Model	Outputs
Non-Invasive Blood Pressure (NIBP)	A&D Medical	UA-767PC	Systolic Pressure (mmHg) Diastolic Pressure (mmHg) Pulse Rate (per minute)
Weight Scale	A&D Medical	UC-321 PL	Weight (up 450 lbs. or 205 kg.)
Glucometer	J&J Lifescan	One Touch Ultra	Blood Glucose (mg/dL)
Extra device port for future expansion			

\* Ask us if you would like support for alternative devices.